


Imaging of ballistic wounds, bullet composition and implications for MRI safety

December 30 2020

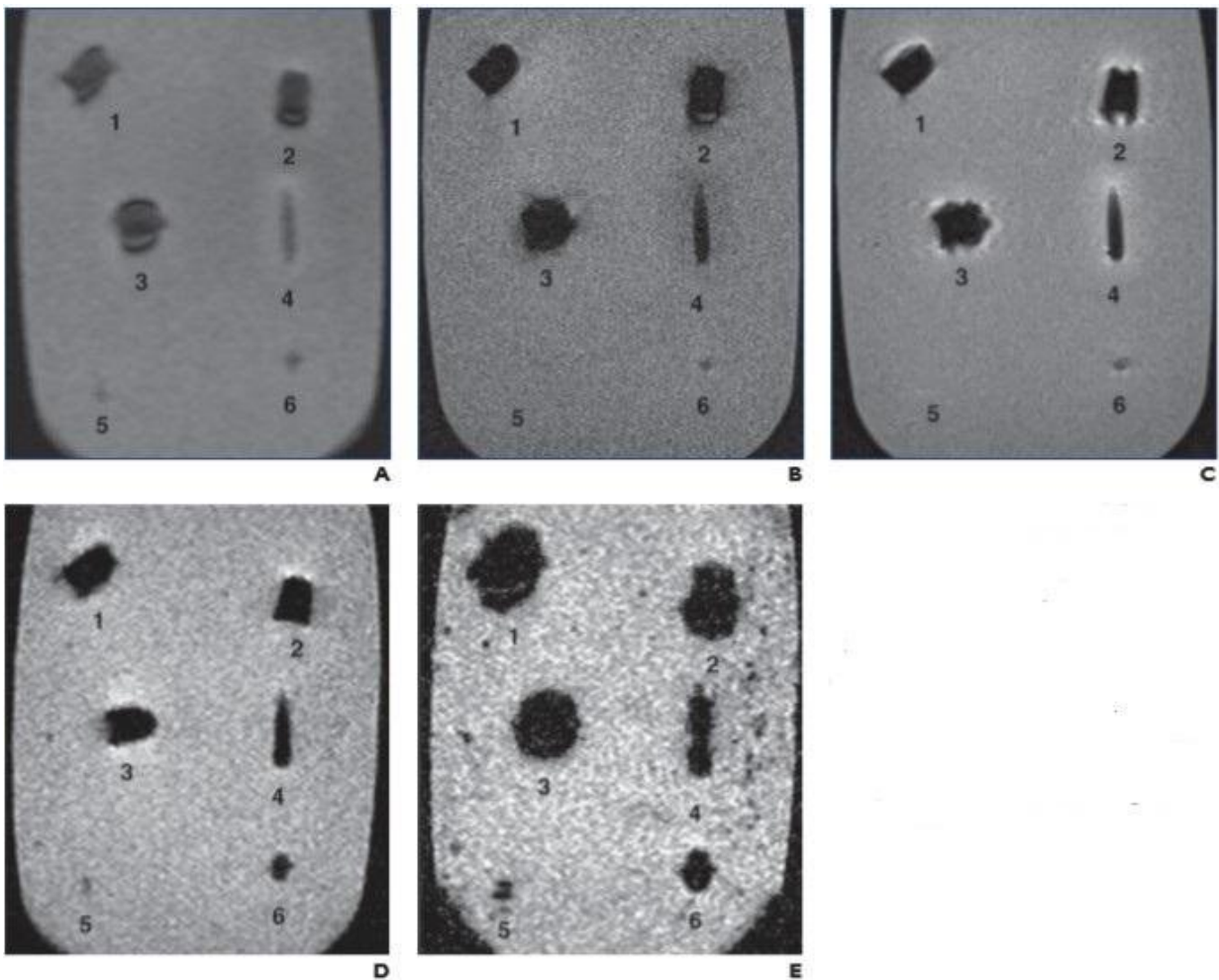

Scout (A), T1-weighted spin-echo (SE) (B), T2-weighted SE (C), T2-weighted gradient-recalled echo (GRE) (TR/TE, 500/10; D), and T2-weighted GRE (TR/TE, 700/30; E) MR images show jacket hollow point .45 automatic Colt pistol bullet (Corbon) (1), solid lead .45 Long Colt bullet (Winchester) (2), full metal jacket (FMJ) automatic Colt pistol bullet (Winchester) (3), 5.56-mm FMJ

bullet (Federal Ammunition) (4), #7 lead shotgun pellet (Winchester) (5), and 5-mm lead air gun pellet (Sheridan) (6). On all sequences, metallic artifact is minimal. Although metallic artifact increases or blooms with increased TR/TE in GRE images (D and E), amount of surrounding distortion is still minimal. Credit: American Roentgen Ray Society (ARRS), *American Journal of Roentgenology* (AJR)


According to an article in ARRS' *American Journal of Roentgenology* (AJR), because patients with ballistic embedded fragments are frequently denied MRI (due to indeterminate bullet composition sans shell casings), radiography and CT can be used to identify nonferromagnetic projectiles that are safe for MRI.

"Commercially available handgun and shotgun ammunition representing projectiles commonly encountered in a [clinical setting](#) was fired into ballistic gelatin as a surrogate for [human tissue](#)," explained first author Arthur J. Fountain from the department of radiology and imaging sciences at Emory University.

After obtaining radiographs and CT images of these gelatin blocks, Fountain and colleagues then obtained MR images of unfired bullets suspended in gelatin blocks using T1- and T2-weighted sequences. Magnetic attractive force, rotational torque, and heating effects of unfired bullets were assessed at 1.5 T.

Based upon debris trail and primary [projectile](#) deformation, the team separated the fired bullets into two groups: ferromagnetic and nonferromagnetic. Although ferromagnetic bullets showed mild torque forces and marked imaging artifacts at 1.5 T, nonferromagnetic bullets did not exhibit these effects.

Importantly, heating above the Food and Drug Administration limit of 2°C was not observed in any of the projectiles tested.

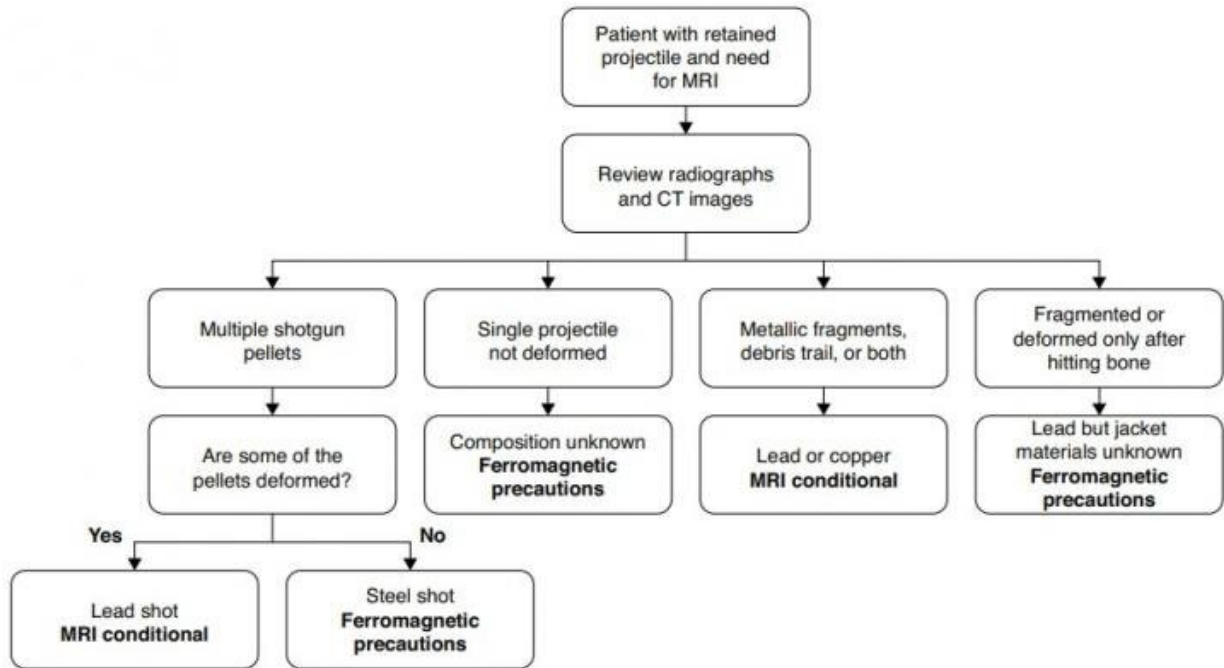

Recommendations for MRI are in bold. MRI conditional indicates imaging is safe at 1.5 T. Ferromagnetic precautions indicate risk-benefit analysis is required before proceeding with MRI. Credit: American Roentgen Ray Society (ARRS), *American Journal of Roentgenology* (AJR)

Additionally, the authors of this AJR article presented a triage algorithm for patients with retained ballistic fragments. "In particular," Fountain et al. described, "a projectile that leaves a metallic debris trail from entry to final position or has been appreciably deformed is of copper, copper-alloy, or lead composition with a partial jacketed configuration or represents lead shotgun shot and does not pose a significant risk for imaging at 1.5 T or less, regardless of when the injury occurred."

"Nonferromagnetic ballistic projectiles do not undergo movement or heating during MRI, and the imaging modality can be performed when medically necessary without undue risk and with limited artifact susceptibility on the resulting [images](#), even when the projectile is in or near a vital structure," the authors concluded.

More information: Arthur J. Fountain et al, Imaging Appearance of Ballistic Wounds Predicts Bullet Composition: Implications for MRI Safety, *American Journal of Roentgenology* (2020). [DOI: 10.2214/AJR.20.23648](#)

Provided by American Roentgen Ray Society

Citation: Imaging of ballistic wounds, bullet composition and implications for MRI safety (2020, December 30) retrieved 16 April 2023 from <https://medicalxpress.com/news/2020-12-imaging-ballistic-wounds-bullet-composition.html>

This document is subject to copyright. Apart from any fair dealing for the purpose of private study or research, no part may be reproduced without the written permission. The content is provided for information purposes only.